
Example Photon Treatment Plan with Direct aperture optimization

Table of Contents

.....	1
Patient Data Import	1
Treatment Plan	1
Generate Beam Geometry STF	2
Dose Calculation	2
Inverse Planning for IMRT	3
Sequencing	7
DAO - Direct Aperture Optimization	8
Aperture visualization	10
Indicator Calculation and display of DVH and QI	12

%%

Copyright 2017 the matRad development team.

This file is part of the matRad project. It is subject to the license terms in the LICENSE file found in the top-level directory of this distribution and at <https://github.com/e0404/matRad/LICENSES.txt>. No part of the matRad project, including this file, may be copied, modified, propagated, or distributed except according to the terms contained in the LICENSE file.

%%

In this example we will show (i) how to load patient data into matRad (ii) how to setup a photon dose calculation and (iii) how to inversely optimize directly from command window in MatLab. (iv) how to apply a sequencing algorithm (v) how to run a direct aperture optimization (iv) how to visually and quantitatively evaluate the result

Patient Data Import

Let's begin with a clear Matlab environment and import the head & neck patient into your workspace.

```
clc,clear,close all;  
load('HEAD_AND_NECK.mat');
```

Treatment Plan

The next step is to define your treatment plan labeled as 'pln'. This structure requires input from the treatment planner and defines the most important cornerstones of your treatment plan.

```
pln.radiationMode = 'photons'; % either photons / protons / carbon  
pln.machine       = 'Generic';  
pln.numOfFractions = 30;  
  
pln.propOpt.bioOptimization = 'none';
```

```
pln.propStf.gantryAngles = [0:72:359];  
pln.propStf.couchAngles = [0 0 0 0 0];  
pln.propStf.bixelWidth = 5;  
pln.propStf.numOfBeams = numel(pln.propStf.gantryAngles);  
pln.propStf.isoCenter = ones(pln.propStf.numOfBeams,1) *  
    matRad_getIsoCenter(cst,ct,0);
```

Enable sequencing and direct aperture optimization (DAO).

```
pln.propOpt.runSequencing = 1;  
pln.propOpt.runDAO = 1;
```

Generate Beam Geometry STF

```
stf = matRad_generateStf(ct,cst,pln);
```

```
matRad: Generating stf struct... Progress: 100.00 %
```

Dose Calculation

Lets generate dosimetric information by pre-computing dose influence matrices for unit beamlet intensities. Having dose influences available allows for subsequent inverse optimization.

```
dij = matRad_calcPhotonDose(ct,stf,pln,cst);
```

```
Warning: ray does not hit patient. Trying to fix afterwards...matRad:  
Photon dose calculation...
```

```
Beam 1 of 5:
```

```
matRad: calculate radiological depth cube...done
```

```
SSD = 919mm
```

```
matRad: Uniform primary photon fluence -> pre-compute kernel  
convolution for SSD = 919 mm ...
```

```
Progress: 100.00 %
```

```
Beam 2 of 5:
```

```
matRad: calculate radiological depth cube...done
```

```
SSD = 943mm
```

```
matRad: Uniform primary photon fluence -> pre-compute kernel  
convolution for SSD = 943 mm ...
```

```
Progress: 100.00 %
```

```
Beam 3 of 5:
```

```
matRad: calculate radiological depth cube...done
```

```
SSD = 928mm
```

```
matRad: Uniform primary photon fluence -> pre-compute kernel  
convolution for SSD = 928 mm ...
```

```
Progress: 100.00 %
```

```
Beam 4 of 5:
```

```
matRad: calculate radiological depth cube...done
```

```
SSD = 908mm
```

```
matRad: Uniform primary photon fluence -> pre-compute kernel  
convolution for SSD = 908 mm ...
```

```
Progress: 100.00 %
```

```
Beam 5 of 5:
```

```
matRad: calculate radiological depth cube...done
```

```
SSD = 934mm
```

```
matRad: Uniform primary photon fluence -> pre-compute kernel
convolution for SSD = 934 mm ...
Progress: 100.00 %
```

Inverse Planning for IMRT

The goal of the fluence optimization is to find a set of beamlet weights which yield the best possible dose distribution according to the predefined clinical objectives and constraints underlying the radiation treatment. Once the optimization has finished, trigger once the GUI to visualize the optimized dose cubes.

```
resultGUI = matRad_fluenceOptimization(dij,cst,pln);
matRadGUI;
```

```
*****
This program contains Ipopt, a library for large-scale nonlinear
optimization.
Ipopt is released as open source code under the Eclipse Public
License (EPL).
For more information visit http://projects.coin-or.org/Ipopt
*****
```

This is Ipopt version 3.12.4, running with linear solver ma57.

```
Number of nonzeros in equality constraint Jacobian...:      0
Number of nonzeros in inequality constraint Jacobian.:      0
Number of nonzeros in Lagrangian Hessian.....:          0
```

```
Total number of variables.....:      5154
      variables with only lower bounds:      5154
      variables with lower and upper bounds:      0
      variables with only upper bounds:      0
Total number of equality constraints.....:      0
Total number of inequality constraints.....:      0
      inequality constraints with only lower bounds:      0
      inequality constraints with lower and upper bounds:      0
      inequality constraints with only upper bounds:      0
```

iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du	alpha_pr	ls
0	5.9484151e+02	0.00e+00	8.83e+00	0.0	0.00e+00	-	0.00e+00	0.00e+00	0
1	2.8033240e+02	0.00e+00	7.59e+00	-1.1	3.36e+00	-	8.61e-01	1.57e-01f	1
2	2.1657816e+02	0.00e+00	2.23e+00	-0.6	1.56e-01	-	9.36e-01	1.00e+00f	1
3	1.8482123e+02	0.00e+00	1.54e+00	-1.7	8.68e-02	-	9.97e-01	1.00e+00f	1
4	1.5093174e+02	0.00e+00	9.50e-01	-2.3	2.19e-01	-	9.98e-01	6.37e-01f	1
5	1.3090358e+02	0.00e+00	1.29e+00	-3.0	2.43e-01	-	1.00e+00	4.02e-01f	1
6	1.1794004e+02	0.00e+00	2.25e+00	-3.4	3.97e-01	-	1.00e+00	2.15e-01f	1

Example Photon Treatment Plan
with Direct aperture optimization

7	1.0099961e+02	0.00e+00	1.17e+00	-2.8	3.80e-01	-	1.00e+00
3.55e-01f	1						
8	9.1986942e+01	0.00e+00	1.05e+00	-3.1	3.67e-01	-	1.00e+00
2.05e-01f	1						
9	8.3502518e+01	0.00e+00	1.51e+00	-2.1	3.05e-01	-	9.85e-01
2.81e-01f	1						
iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du
alpha_pr	ls						
10	4.9306911e+02	0.00e+00	7.73e+00	-0.1	1.14e+01	-	1.48e-01
1.53e-01f	1						
11	1.4619937e+02	0.00e+00	2.96e+00	-0.9	9.84e-01	-	1.00e+00
1.00e+00f	1						
12	1.3768423e+02	0.00e+00	1.32e+00	-0.9	1.00e-01	-	1.00e+00
1.00e+00f	1						
13	9.0840902e+01	0.00e+00	4.47e-01	-1.5	3.28e-01	-	8.13e-01
1.00e+00f	1						
14	7.7213074e+01	0.00e+00	6.28e-01	-1.5	1.74e-01	-	1.00e+00
1.00e+00f	1						
15	6.7525359e+01	0.00e+00	4.29e-01	-2.3	2.07e-01	-	9.97e-01
7.91e-01f	1						
16	6.5610126e+01	0.00e+00	1.39e+00	-3.4	9.70e-02	-	1.00e+00
2.74e-01f	1						
17	6.0607242e+01	0.00e+00	5.24e-01	-3.6	1.75e-01	-	9.99e-01
6.09e-01f	1						
18	5.8552703e+01	0.00e+00	3.16e-01	-2.8	1.20e-01	-	1.00e+00
4.00e-01f	1						
19	5.7455023e+01	0.00e+00	1.01e+00	-3.2	8.25e-02	-	1.00e+00
3.97e-01f	1						
iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du
alpha_pr	ls						
20	5.6202297e+01	0.00e+00	3.26e-01	-2.4	5.95e-02	-	9.54e-01
1.00e+00f	1						
21	5.5412369e+01	0.00e+00	3.79e-01	-2.4	7.71e-02	-	7.10e-01
9.17e-01f	1						
22	5.5007773e+01	0.00e+00	1.06e+00	-3.3	5.63e-02	-	1.00e+00
2.05e-01f	1						
23	5.4174549e+01	0.00e+00	6.75e-01	-9.2	8.73e-02	-	8.09e-01
2.90e-01f	1						
24	5.3351574e+01	0.00e+00	4.86e-01	-3.0	6.42e-02	-	6.34e-01
4.67e-01f	1						
25	5.3138334e+01	0.00e+00	4.08e-01	-2.6	5.20e-02	-	5.58e-01
1.00e+00f	1						
26	5.2671186e+01	0.00e+00	2.44e-01	-2.8	3.44e-02	-	6.66e-01
7.04e-01f	1						
27	5.2125268e+01	0.00e+00	4.11e-01	-3.0	7.38e-02	-	9.98e-01
3.37e-01f	1						
28	5.1152934e+01	0.00e+00	2.91e-01	-3.5	1.50e-01	-	1.00e+00
3.27e-01f	1						
29	5.0837115e+01	0.00e+00	2.76e-01	-2.8	7.10e-02	-	6.61e-01
2.89e-01f	1						
iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du
alpha_pr	ls						
30	5.0171847e+01	0.00e+00	3.64e-01	-3.2	1.12e-01	-	8.17e-01
3.77e-01f	1						

Example Photon Treatment Plan
with Direct aperture optimization

31	4.9820321e+01	0.00e+00	3.59e-01	-3.1	8.04e-02	-	6.70e-01
	3.28e-01f	1					
32	4.9426185e+01	0.00e+00	2.82e-01	-9.1	1.32e-01	-	5.66e-01
	2.18e-01f	1					
33	4.9223616e+01	0.00e+00	3.13e-01	-4.5	1.23e-01	-	7.97e-01
	1.22e-01f	1					
34	4.8700424e+01	0.00e+00	2.18e-01	-4.9	1.62e-01	-	1.00e+00
	2.68e-01f	1					
35	4.8478784e+01	0.00e+00	4.06e-01	-3.8	7.83e-02	-	1.00e+00
	2.34e-01f	1					
36	4.8177447e+01	0.00e+00	2.95e-01	-4.0	1.02e-01	-	8.57e-01
	2.80e-01f	1					
37	4.8001419e+01	0.00e+00	3.51e-01	-4.5	1.22e-01	-	6.51e-01
	1.50e-01f	1					
38	4.7943712e+01	0.00e+00	6.39e-01	-6.1	8.55e-02	-	5.68e-01
	6.86e-02f	1					
39	4.7713931e+01	0.00e+00	3.21e-01	-4.6	1.25e-01	-	6.92e-01
	2.10e-01f	1					
iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du
	alpha_pr	ls					
40	4.7502222e+01	0.00e+00	3.51e-01	-4.9	1.21e-01	-	8.48e-01
	2.18e-01f	1					
41	4.7407062e+01	0.00e+00	3.43e-01	-10.4	8.33e-02	-	4.94e-01
	1.58e-01f	1					
42	4.7257084e+01	0.00e+00	2.10e-01	-5.5	1.55e-01	-	9.82e-01
	1.57e-01f	1					
43	4.7185505e+01	0.00e+00	5.08e-01	-5.7	9.98e-02	-	7.65e-01
	1.15e-01f	1					
44	4.7032088e+01	0.00e+00	2.44e-01	-4.9	2.11e-01	-	4.84e-01
	1.47e-01f	1					
45	4.7018421e+01	0.00e+00	6.68e-01	-10.6	7.75e-02	-	5.81e-01
	3.08e-02f	1					
46	4.6910455e+01	0.00e+00	3.10e-01	-5.4	1.47e-01	-	7.37e-01
	1.35e-01f	1					
47	4.6833803e+01	0.00e+00	2.91e-01	-5.6	1.36e-01	-	4.73e-01
	1.07e-01f	1					
48	4.6787414e+01	0.00e+00	5.59e-01	-10.8	9.09e-02	-	5.80e-01
	9.84e-02f	1					
49	4.6724826e+01	0.00e+00	5.12e-01	-4.7	8.47e-02	-	6.70e-01
	1.54e-01f	1					
iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du
	alpha_pr	ls					
50	4.6661660e+01	0.00e+00	2.85e-01	-4.1	1.07e-01	-	4.50e-01
	1.45e-01f	1					
51	4.6640277e+01	0.00e+00	1.98e-01	-10.3	1.36e-01	-	3.91e-01
	3.79e-02f	1					
52	4.6528248e+01	0.00e+00	1.82e-01	-4.5	1.85e-01	-	9.18e-01
	1.53e-01f	1					
53	4.6513717e+01	0.00e+00	6.31e-01	-4.2	8.69e-02	-	7.12e-01
	4.58e-02f	1					
54	4.6441858e+01	0.00e+00	2.45e-01	-5.1	1.25e-01	-	5.68e-01
	1.53e-01f	1					
55	4.6462436e+01	0.00e+00	1.65e-01	-3.4	1.31e-01	-	2.82e-01
	2.49e-01f	1					

Example Photon Treatment Plan
with Direct aperture optimization

```

56  4.6454990e+01 0.00e+00 4.17e-01 -3.9 3.25e-02 - 6.61e-01
5.59e-02f 1
57  4.6409093e+01 0.00e+00 3.35e-01 -3.9 6.57e-02 - 4.65e-01
1.80e-01f 1

```

Number of Iterations.....: 57

	(scaled)	(unscaled)
Objective.....:	4.6409092889633186e+01	
	4.6409092889633186e+01	
Dual infeasibility.....:	3.3452863040260739e-01	
	3.3452863040260739e-01	
Constraint violation.....:	0.0000000000000000e+00	
	0.0000000000000000e+00	
Complementarity.....:	3.8091028351478075e-04	
	3.8091028351478075e-04	
Overall NLP error.....:	3.3452863040260739e-01	
	3.3452863040260739e-01	

Number of objective function evaluations	= 58
Number of objective gradient evaluations	= 58
Number of equality constraint evaluations	= 0
Number of inequality constraint evaluations	= 0
Number of equality constraint Jacobian evaluations	= 0
Number of inequality constraint Jacobian evaluations	= 0
Number of Lagrangian Hessian evaluations	= 0
Total CPU secs in IPOPT (w/o function evaluations)	= 1.924
Total CPU secs in NLP function evaluations	= 18.540

EXIT: Solved To Acceptable Level.



*** IPOPT DONE ***

Calculating final cubes...

Warning: 'popupmenu' control requires that 'Value' be an integer
within

Character vector range

Control will not be rendered until all of its parameter values are
valid

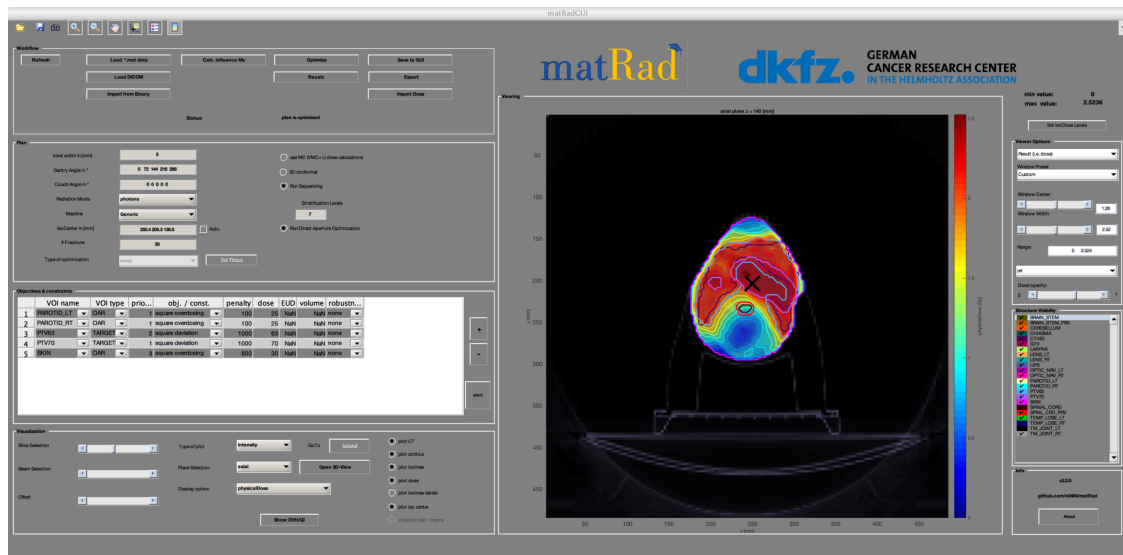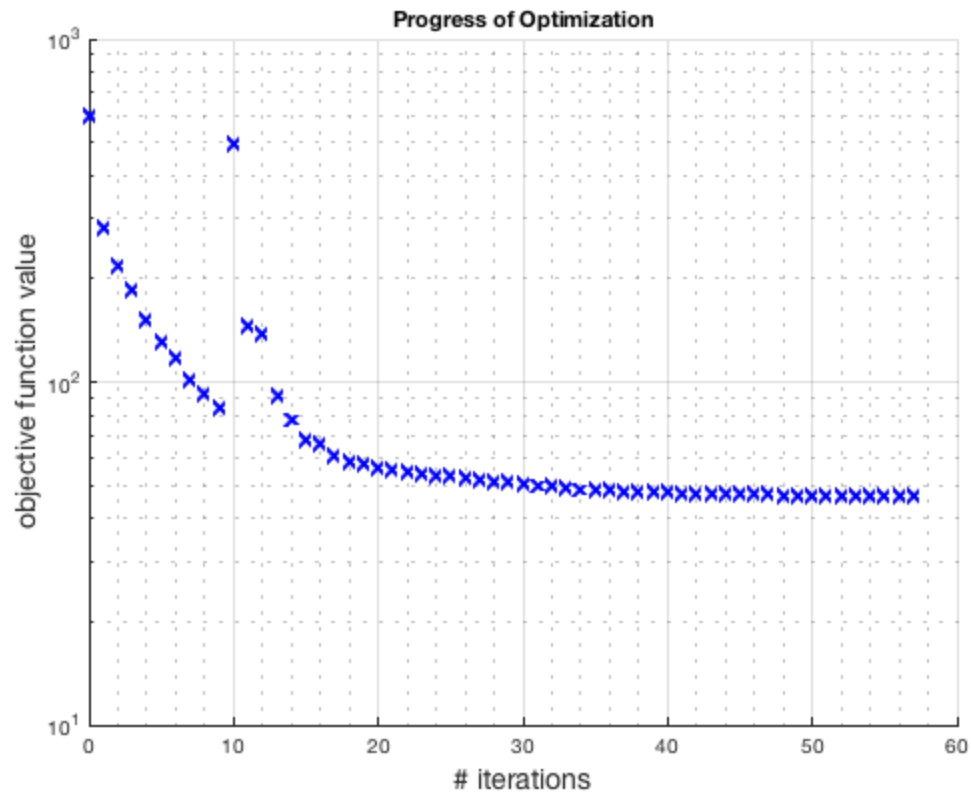

Sequencing

This is a multileaf collimator leaf sequencing algorithm that is used in order to modulate the intensity of the beams with multiple static segments, so that translates each intensity map into a set of deliverable aperture shapes.

```
resultGUI = matRad_siochiLeafSequencing(resultGUI,stf,dij,5);
```

DAO - Direct Aperture Optimization

The Direct Aperture Optimization is an optimization approach where we directly optimize aperture shapes and weights.

```
resultGUI =
  matRad_directApertureOptimization(dij,cst,resultGUI.apertureInfo,resultGUI,pln);
```

```
*****
This program contains Ipopt, a library for large-scale nonlinear
optimization.
Ipopt is released as open source code under the Eclipse Public
License (EPL).
For more information visit http://projects.coin-or.org/Ipopt
*****
```

This is Ipopt version 3.12.4, running with linear solver ma57.

```
Number of nonzeros in equality constraint Jacobian...:      0
Number of nonzeros in inequality constraint Jacobian.:    7304
Number of nonzeros in Lagrangian Hessian.....:          0
```

```
Total number of variables.....:    7393
      variables with only lower bounds:      89
      variables with lower and upper bounds:  7304
      variables with only upper bounds:       0
Total number of equality constraints.....:      0
Total number of inequality constraints.....:    3652
      inequality constraints with only lower bounds:  3652
      inequality constraints with lower and upper bounds:  0
      inequality constraints with only upper bounds:  0
```

iter	objective	inf_pr	inf_du	lg(mu)	d	lg(rg)	alpha_du	alpha_pr	ls
0	5.5447951e+01	0.00e+00	2.45e+01	0.0	0.00e+00	-	0.00e+00	0.00e+00	0
1	2.6135015e+05	0.00e+00	1.44e+04	0.8	1.15e+01	-	1.00e+00	3.68e-01h	1
2	3.0441990e+02	0.00e+00	3.36e+02	0.7	3.92e+00	-	1.00e+00	1.00e+00f	1
3	1.8933618e+02	0.00e+00	1.79e+02	-0.6	5.01e-02	-	1.00e+00	1.00e+00f	1
4	9.3170671e+01	0.00e+00	6.49e+01	-1.3	7.80e-02	-	1.00e+00	1.00e+00f	1
5	7.2754136e+01	0.00e+00	3.58e+01	-2.0	3.25e-02	-	9.99e-01	1.00e+00f	1
6	6.1071647e+01	0.00e+00	3.15e+01	-3.1	3.95e-02	-	1.00e+00	1.00e+00f	1
7	5.9424773e+01	0.00e+00	6.40e+01	-3.7	4.26e-02	-	1.00e+00	1.00e+00f	1
8	5.6404373e+01	0.00e+00	1.46e+01	-4.8	1.43e-02	-	1.00e+00	1.00e+00f	1

Example Photon Treatment Plan
with Direct aperture optimization

```

    9  5.6023161e+01  0.00e+00  7.77e+00  -6.4  4.42e-03  -  1.00e+00
1.00e+00f  1
iter    objective    inf_pr    inf_du lg(mu)  ||d||  lg(rg) alpha_du
alpha_pr  ls
   10  5.5505052e+01  0.00e+00  1.06e+01  -7.9  9.10e-03  -  1.00e+00
1.00e+00f  1
   11  5.5017729e+01  0.00e+00  7.54e+00  -9.4  1.11e-02  -  1.00e+00
1.00e+00f  1
   12  5.4722128e+01  0.00e+00  1.59e+01 -10.7  2.50e-02  -  1.00e+00
1.00e+00f  1
   13  5.4468243e+01  0.00e+00  1.01e+01 -11.0  1.35e-02  -  1.00e+00
1.00e+00f  1
   14  5.4266782e+01  0.00e+00  3.02e+00 -11.0  8.48e-03  -  1.00e+00
1.00e+00f  1
   15  5.4208602e+01  0.00e+00  2.60e+00 -11.0  2.11e-03  -  1.00e+00
1.00e+00f  1
   16  5.4026245e+01  0.00e+00  4.51e+00 -11.0  1.22e-02  -  1.00e+00
1.00e+00f  1
   17  5.3974784e+01  0.00e+00  5.74e+00 -11.0  1.85e-02  -  1.00e+00
2.50e-01f  3
   18  5.3916667e+01  0.00e+00  3.73e+00 -11.0  6.07e-03  -  1.00e+00
1.00e+00f  1
   19  5.3841410e+01  0.00e+00  3.38e+00 -11.0  8.42e-03  -  1.00e+00
1.00e+00f  1
iter    objective    inf_pr    inf_du lg(mu)  ||d||  lg(rg) alpha_du
alpha_pr  ls
   20  5.3778191e+01  0.00e+00  3.69e+00 -11.0  6.97e-03  -  1.00e+00
1.00e+00f  1
   21  5.3706229e+01  0.00e+00  1.64e+00 -11.0  8.95e-03  -  1.00e+00
1.00e+00f  1
   22  5.3693382e+01  0.00e+00  4.90e+00 -11.0  7.71e-03  -  1.00e+00
1.00e+00f  1
   23  5.3654180e+01  0.00e+00  1.56e+00 -11.0  3.59e-03  -  1.00e+00
1.00e+00f  1
   24  5.3643953e+01  0.00e+00  1.61e+00 -11.0  1.52e-03  -  1.00e+00
1.00e+00f  1

```

Number of Iterations.....: 24


	(scaled)	(unscaled)
Objective.....:	5.3643952996389672e+01	
	5.3643952996389672e+01	
Dual infeasibility.....:	1.6101521274568646e+00	
	1.6101521274568646e+00	
Constraint violation.....:	0.0000000000000000e+00	
	0.0000000000000000e+00	
Complementarity.....:	1.0000000001370447e-11	
	1.0000000001370447e-11	
Overall NLP error.....:	1.6101521274568646e+00	
	1.6101521274568646e+00	

Number of objective function evaluations = 31
Number of objective gradient evaluations = 25

Number of equality constraint evaluations	=	0
Number of inequality constraint evaluations	=	31
Number of equality constraint Jacobian evaluations	=	0
Number of inequality constraint Jacobian evaluations	=	25
Number of Lagrangian Hessian evaluations	=	0
Total CPU secs in IPOPT (w/o function evaluations)	=	1.423
Total CPU secs in NLP function evaluations	=	14.050

EXIT: Solved To Acceptable Level.

*** IPOPT DONE ***

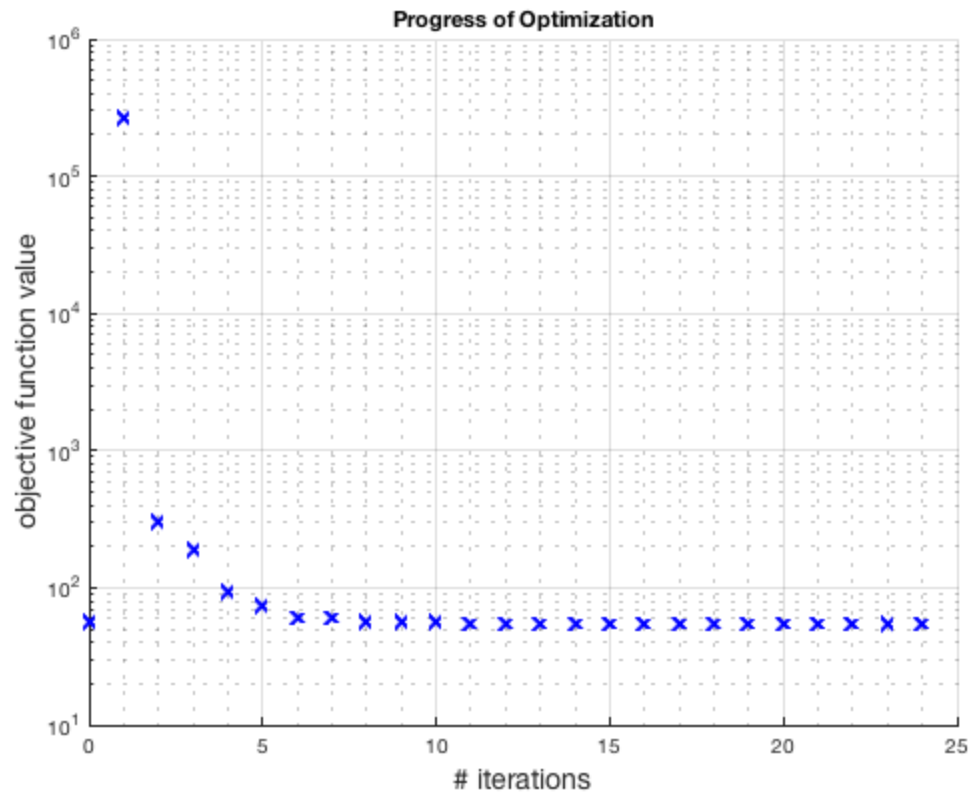



Aperture visualization

Use a matrad function to visualize the resulting aperture shapes

```
matRad_visApertureInfo(resultGUI.apertureInfo);
```

Example Photon Treatment Plan with Direct aperture optimization

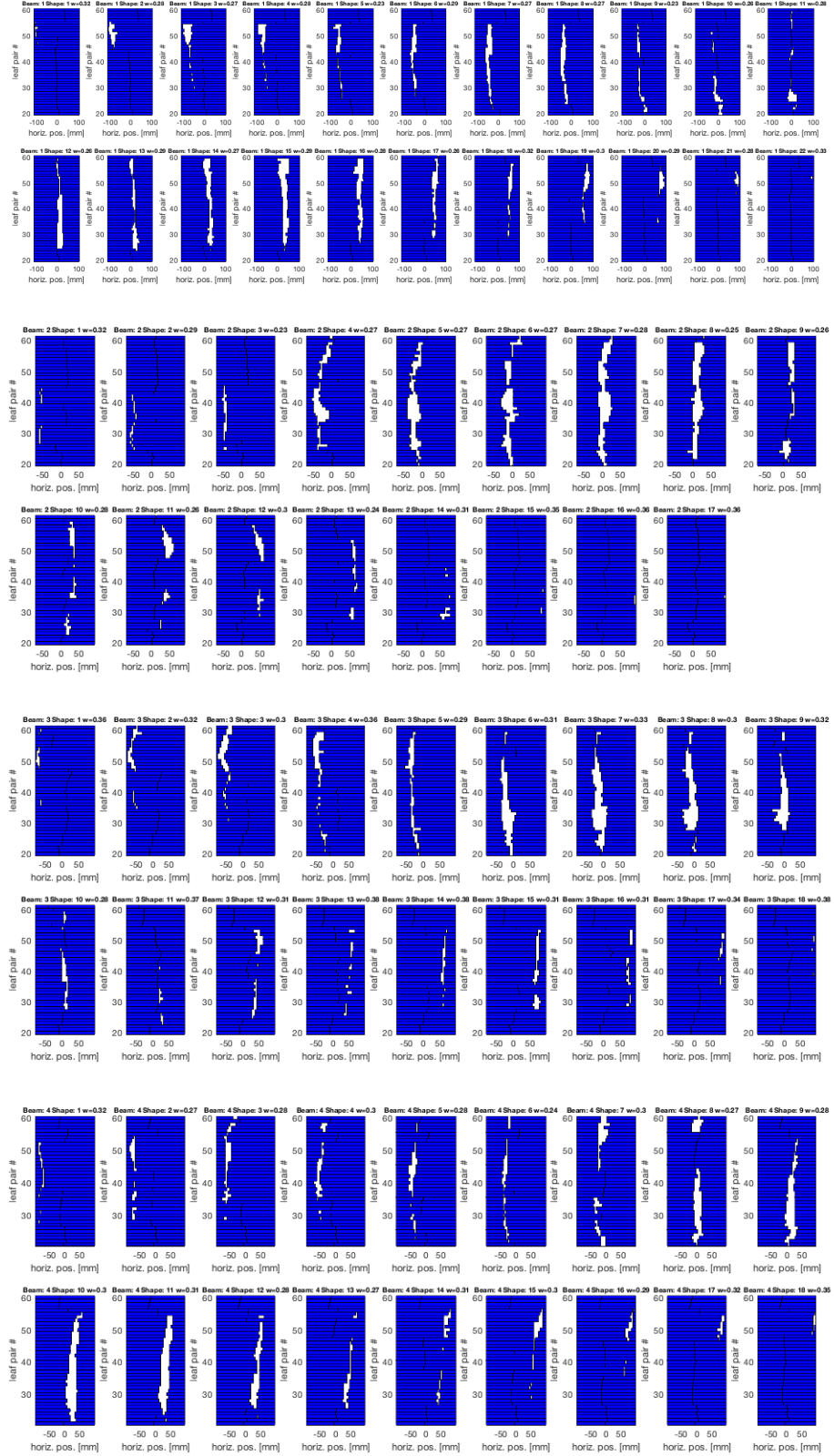

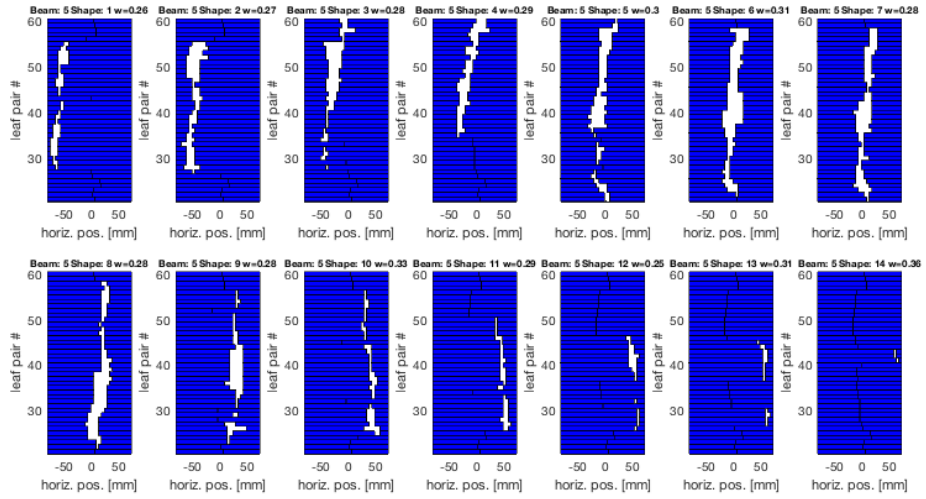

Indicator Calculation and display of DVH and QI

```
[dvh,qi] = matRad_indicatorWrapper(cst,pln,resultGUI);
```

```
0          BRAIN_STEM - Mean dose = 0.37 Gy +/- 0.31 Gy (Max dose
= 1.24 Gy, Min dose = 0.03 Gy)
          D2% = 1.06 Gy, D5% = 0.95 Gy, D50% =
0.29 Gy, D95% = 0.05 Gy, D98% = 0.04 Gy,
          V0Gy = 100.00%, V0.5Gy = 29.83%, V1Gy =
3.58%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,
```

```
1          BRAIN_STEM_PRV - Mean dose = 0.37 Gy +/- 0.33 Gy (Max dose
= 1.27 Gy, Min dose = 0.02 Gy)
          D2% = 1.08 Gy, D5% = 0.98 Gy, D50% =
0.29 Gy, D95% = 0.03 Gy, D98% = 0.03 Gy,
          V0Gy = 100.00%, V0.5Gy = 31.48%, V1Gy =
4.15%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,
```

```
2          CEREBELLUM - Mean dose = 0.59 Gy +/- 0.31 Gy (Max dose
= 1.66 Gy, Min dose = 0.03 Gy)
          D2% = 1.32 Gy, D5% = 1.13 Gy, D50% =
0.58 Gy, D95% = 0.11 Gy, D98% = 0.06 Gy,
          V0Gy = 100.00%, V0.5Gy = 58.88%, V1Gy =
10.88%, V1.5Gy = 0.29%, V2Gy = 0.00%, V2.5Gy = 0.00%,
```

```
3          CHIASMA - Mean dose = 0.10 Gy +/- 0.03 Gy (Max dose
= 0.13 Gy, Min dose = 0.06 Gy)
          D2% = 0.13 Gy, D5% = 0.13 Gy, D50% =
0.11 Gy, D95% = 0.06 Gy, D98% = 0.06 Gy,
          V0Gy = 100.00%, V0.5Gy = 0.00%, V1Gy =
0.00%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,
```

5 CTV63 - Mean dose = 2.14 Gy +/- 0.15 Gy (Max dose = 2.53 Gy, Min dose = 0.90 Gy)
D2% = 2.42 Gy, D5% = 2.37 Gy, D50% = 2.14 Gy, D95% = 1.92 Gy, D98% = 1.78 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy = 99.98%, V1.5Gy = 99.56%, V2Gy = 88.83%, V2.5Gy = 0.14%,
Warning: target has no objective that penalizes underdosage,

6 GTV - Mean dose = 2.34 Gy +/- 0.06 Gy (Max dose = 2.57 Gy, Min dose = 2.13 Gy)
D2% = 2.47 Gy, D5% = 2.45 Gy, D50% = 2.34 Gy, D95% = 2.24 Gy, D98% = 2.21 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy = 100.00%, V1.5Gy = 100.00%, V2Gy = 100.00%, V2.5Gy = 0.52%,
Warning: target has no objective that penalizes underdosage,

7 LARYNX - Mean dose = 1.24 Gy +/- 0.23 Gy (Max dose = 1.85 Gy, Min dose = 0.89 Gy)
D2% = 1.75 Gy, D5% = 1.67 Gy, D50% = 1.19 Gy, D95% = 0.95 Gy, D98% = 0.92 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy = 85.03%, V1.5Gy = 15.65%, V2Gy = 0.00%, V2.5Gy = 0.00%,

8 LENS_LT - Mean dose = 0.01 Gy +/- 0.00 Gy (Max dose = 0.01 Gy, Min dose = 0.00 Gy)
D2% = 0.01 Gy, D5% = 0.01 Gy, D50% = 0.01 Gy, D95% = 0.00 Gy, D98% = 0.00 Gy,
V0Gy = 100.00%, V0.5Gy = 0.00%, V1Gy = 0.00%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

9 LENS_RT - Mean dose = 0.01 Gy +/- 0.00 Gy (Max dose = 0.01 Gy, Min dose = 0.00 Gy)
D2% = 0.01 Gy, D5% = 0.01 Gy, D50% = 0.01 Gy, D95% = 0.01 Gy, D98% = 0.00 Gy,
V0Gy = 100.00%, V0.5Gy = 0.00%, V1Gy = 0.00%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

10 LIPS - Mean dose = 0.97 Gy +/- 0.38 Gy (Max dose = 2.04 Gy, Min dose = 0.25 Gy)
D2% = 1.75 Gy, D5% = 1.62 Gy, D50% = 0.97 Gy, D95% = 0.40 Gy, D98% = 0.34 Gy,
V0Gy = 100.00%, V0.5Gy = 89.78%, V1Gy = 45.70%, V1.5Gy = 10.75%, V2Gy = 1.08%, V2.5Gy = 0.00%,

11 OPTIC_NRV_LT - Mean dose = 0.06 Gy +/- 0.03 Gy (Max dose = 0.14 Gy, Min dose = 0.02 Gy)
D2% = 0.14 Gy, D5% = 0.14 Gy, D50% = 0.05 Gy, D95% = 0.02 Gy, D98% = 0.02 Gy,
V0Gy = 100.00%, V0.5Gy = 0.00%, V1Gy = 0.00%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

12 OPTIC_NRV_RT - Mean dose = 0.05 Gy +/- 0.03 Gy (Max dose = 0.12 Gy, Min dose = 0.01 Gy)

Example Photon Treatment Plan
with Direct aperture optimization

D2% = 0.12 Gy, D5% = 0.11 Gy, D50% =
0.04 Gy, D95% = 0.02 Gy, D98% = 0.01 Gy,
V0Gy = 100.00%, V0.5Gy = 0.00%, V1Gy =
0.00%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

13 PAROTID_LT - Mean dose = 0.66 Gy +/- 0.28 Gy (Max dose
= 1.65 Gy, Min dose = 0.23 Gy)
D2% = 1.34 Gy, D5% = 1.21 Gy, D50% =
0.62 Gy, D95% = 0.30 Gy, D98% = 0.27 Gy,
V0Gy = 100.00%, V0.5Gy = 68.02%, V1Gy =
12.16%, V1.5Gy = 0.45%, V2Gy = 0.00%, V2.5Gy = 0.00%,

14 PAROTID_RT - Mean dose = 0.67 Gy +/- 0.25 Gy (Max dose
= 1.56 Gy, Min dose = 0.29 Gy)
D2% = 1.28 Gy, D5% = 1.17 Gy, D50% =
0.61 Gy, D95% = 0.37 Gy, D98% = 0.33 Gy,
V0Gy = 100.00%, V0.5Gy = 67.97%, V1Gy =
11.15%, V1.5Gy = 0.16%, V2Gy = 0.00%, V2.5Gy = 0.00%,

16 PTV63 - Mean dose = 2.11 Gy +/- 0.19 Gy (Max dose
= 2.55 Gy, Min dose = 0.62 Gy)
D2% = 2.42 Gy, D5% = 2.37 Gy, D50% =
2.12 Gy, D95% = 1.83 Gy, D98% = 1.70 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy =
99.68%, V1.5Gy = 98.97%, V2Gy = 79.70%, V2.5Gy = 0.13%,
CI = 0.7350, HI = 25.77 for reference dose
of 2.1 Gy

17 PTV70 - Mean dose = 2.32 Gy +/- 0.09 Gy (Max dose
= 2.57 Gy, Min dose = 1.49 Gy)
D2% = 2.47 Gy, D5% = 2.44 Gy, D50% =
2.32 Gy, D95% = 2.18 Gy, D98% = 2.13 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy =
100.00%, V1.5Gy = 99.98%, V2Gy = 99.46%, V2.5Gy = 0.58%,
CI = 0.6456, HI = 11.16 for reference dose
of 2.3 Gy

18 SKIN - Mean dose = 0.54 Gy +/- 0.70 Gy (Max dose
= 2.57 Gy, Min dose = 0.00 Gy)
D2% = 2.28 Gy, D5% = 2.13 Gy, D50% =
0.16 Gy, D95% = 0.00 Gy, D98% = 0.00 Gy,
V0Gy = 100.00%, V0.5Gy = 35.40%, V1Gy =
22.03%, V1.5Gy = 13.25%, V2Gy = 8.05%, V2.5Gy = 0.01%,

19 SPINAL_CORD - Mean dose = 0.79 Gy +/- 0.18 Gy (Max dose
= 1.30 Gy, Min dose = 0.37 Gy)
D2% = 1.16 Gy, D5% = 1.10 Gy, D50% =
0.78 Gy, D95% = 0.55 Gy, D98% = 0.41 Gy,
V0Gy = 100.00%, V0.5Gy = 96.00%, V1Gy =
13.88%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

20 SPINL_CRD_PRV - Mean dose = 0.78 Gy +/- 0.26 Gy (Max dose
= 1.48 Gy, Min dose = 0.07 Gy)

Example Photon Treatment Plan
with Direct aperture optimization

D2% = 1.30 Gy, D5% = 1.20 Gy, D50% =
0.77 Gy, D95% = 0.31 Gy, D98% = 0.21 Gy,
V0Gy = 100.00%, V0.5Gy = 85.47%, V1Gy =
19.80%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,


21 TEMP_LOBE_LT - Mean dose = 0.13 Gy +/- 0.19 Gy (Max dose
= 1.70 Gy, Min dose = 0.01 Gy)
D2% = 0.82 Gy, D5% = 0.49 Gy, D50% =
0.06 Gy, D95% = 0.02 Gy, D98% = 0.01 Gy,
V0Gy = 100.00%, V0.5Gy = 4.83%, V1Gy =
0.88%, V1.5Gy = 0.10%, V2Gy = 0.00%, V2.5Gy = 0.00%,

22 TEMP_LOBE_RT - Mean dose = 0.10 Gy +/- 0.13 Gy (Max dose
= 0.75 Gy, Min dose = 0.01 Gy)
D2% = 0.51 Gy, D5% = 0.41 Gy, D50% =
0.05 Gy, D95% = 0.01 Gy, D98% = 0.01 Gy,
V0Gy = 100.00%, V0.5Gy = 2.09%, V1Gy =
0.00%, V1.5Gy = 0.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

23 TM_JOINT_LT - Mean dose = 1.39 Gy +/- 0.35 Gy (Max dose
= 1.99 Gy, Min dose = 0.91 Gy)
D2% = 1.96 Gy, D5% = 1.90 Gy, D50% =
1.37 Gy, D95% = 0.92 Gy, D98% = 0.91 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy =
80.00%, V1.5Gy = 30.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

24 TM_JOINT_RT - Mean dose = 1.01 Gy +/- 0.41 Gy (Max dose
= 1.60 Gy, Min dose = 0.61 Gy)
D2% = 1.57 Gy, D5% = 1.53 Gy, D50% =
0.89 Gy, D95% = 0.62 Gy, D98% = 0.61 Gy,
V0Gy = 100.00%, V0.5Gy = 100.00%, V1Gy =
40.00%, V1.5Gy = 20.00%, V2Gy = 0.00%, V2.5Gy = 0.00%,

Example Photon Treatment Plan
with Direct aperture optimization

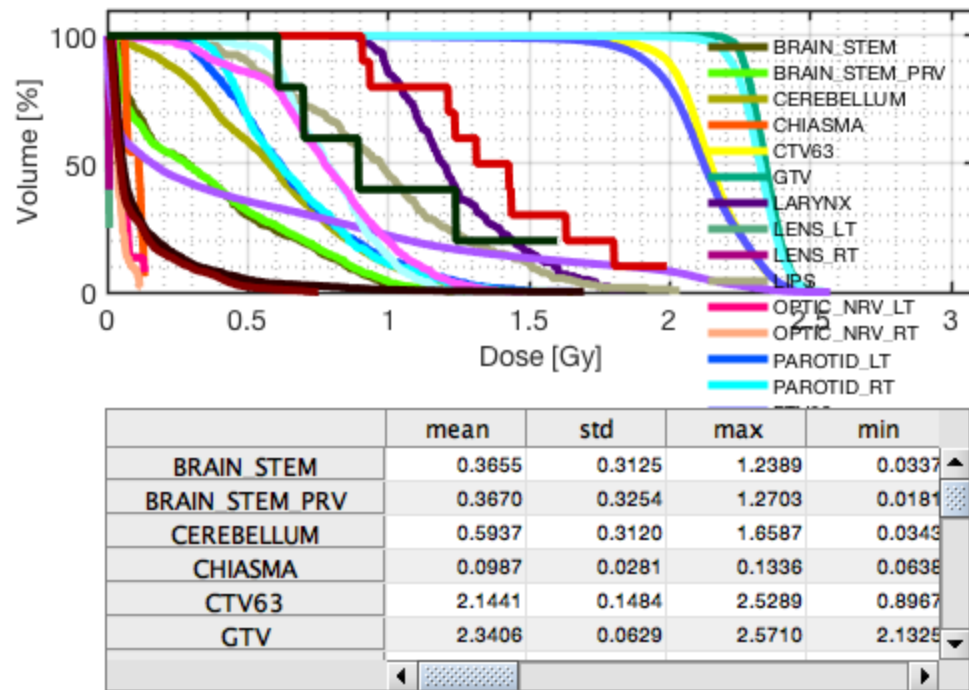

Published with MATLAB® R2017b